

## Refrigerator & Freezer Accessory

The Floor and Wall Bracket Kit is designed to secure your refrigerator\* or freezer to the building floor and/or wall. Helmer refrigerators and freezers include factory threaded mounting holes to simplify installation. Heavy gauge, powder coated steel floor mounting brackets and cabinet anchor hardware are provided. (Bracket Kit is not seismic approved)




### Installation Notes

Kit includes: eyebolts (2), hex nuts (2), brackets (2), hex screws (2)


The hex nuts and eyebolts are for the top holes. The brackets and hex screws are for the bottom holes. Customer to provide floor and/or wall anchor hardware and on-site assembly. Methods for floor/wall anchoring will vary by facility.

### Specifications

Single door refrigerator/freezer (back view)

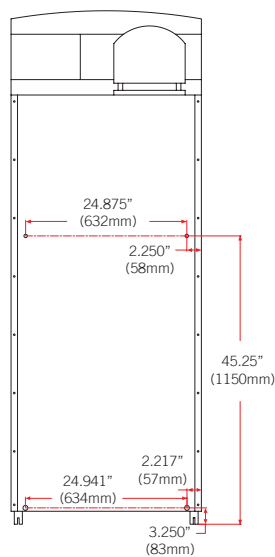

Double door refrigerator/freezer (back view)

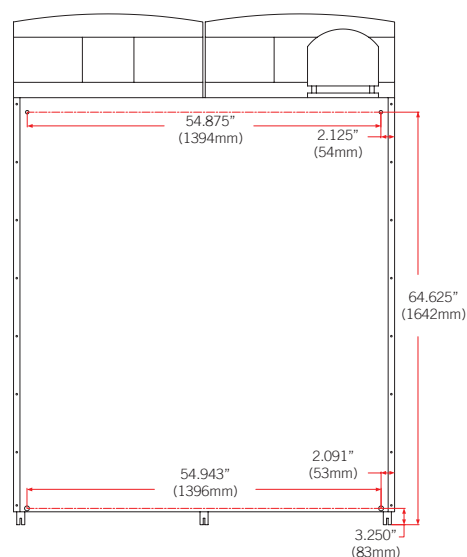

\*Pass-Thru refrigerators include Floor Bracket only.