


# Remote Alarms

For remote location monitoring

Helmer alarms extend your existing system to remote locations in your facility. They are designed for applications where it may be difficult to hear the existing device alarm.

Two models are available. The RA1 offers monitoring of one device or zone and the RA4 provides monitoring of up to four devices or zones.

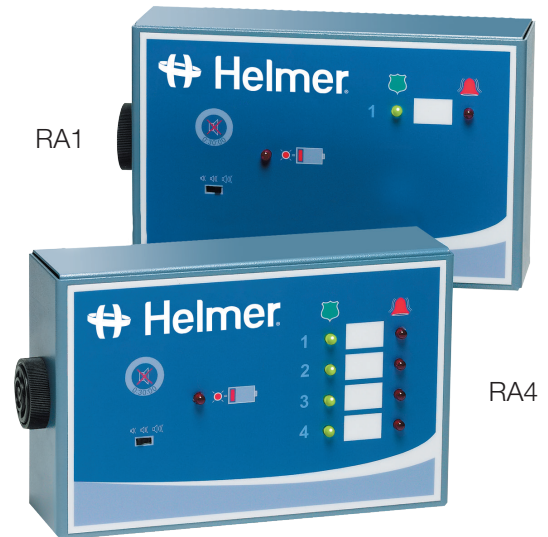

## Features

- Operational at distances up to 3000 feet
- Audible and visual signals indicate an alarm condition
- Countertop, wall or cabinet mountable
- Alarm can be silenced for 30 minutes with audible ring-back
- Compatible with devices featuring dry contact connections
- Adjustable audible alarm volume
- AC powered with non-rechargeable battery back-up
- Low battery indicator
- Each zone independently monitors a normally closed or normally open loop (user selectable)
- Can be connected to a central alarm system
- Certified to applicable UL and CSA standards by a NRTL

## Specifications

Models	RA1	RA4
Capacity	One Device or Zone	Four Devices or Zones
Dimensions (w x h x d)	8.25" x 5.75" x 2.75" 210mm x 146mm x 70mm	
Weight	3 lb / 1.4 kg	
Electrical	100-240V 50/60Hz*	
Wiring	Two-conductor, 22-gauge (not supplied)	
Warranty	2 years parts / 1 year labor	

\* Universal plugs included for USA/Canada and International use.

Specifications subject to change