


Model

i.Series® iLR104-ADA Undercounter Laboratory/Pharmacy Refrigerator

Specifications

Application, Rating and Electrical Data		Controller	
Application	Storage of medical and scientific products	Interface	i.C ³ ® color touchscreen
Storage Volume	4 cu ft (113 L)	Power Switch	On/Off
Temperature Range	2° to 10°C	Controller Type	Microprocessor with touchscreen input and display. Includes USB port for data retrieval.
Electrical Power	115V 60Hz	Security	Lockable door, password protected settings, optional PIN access
Maximum Current	5 FLA	Control Sensor	RTD
Building Supply Rating	15 amp dedicated ground circuit	Communication Ports	(2) USB Ports (1) RJ45 Hub (1) RS232 Port (optional)
Power Plug / Power Cord Length	NEMA 5-15 hospital grade plug (115V) 8-10 ft (2439-3048 mm)	High / Low Alarms	Fully adjustable
Certification / Agency Listing	QPS (Certified to UL and CSA Standards)	High / Low Alarm Test	Automatic (Peltier)
Indoor/Outdoor Use	Indoor use only	Door Ajar Alarm	Fully adjustable
Application Environment	Non-corrosive, non-flammable, non-explosive, 15°C to 32°C (59°F to 90°F)	Power Failure Alarm	Yes
Refrigeration System		Condenser Temperature Alarm	Fully adjustable
Refrigeration System	Forced-air circulation	Event Log	All alarms, door openings, download via USB
Compressor / Number	.17 HP (115V) hermetic / 1	Temperature Log	60 days, download via USB
Condenser Type / Number	Air-cooled fin and tube / 1	Access Log	60 days, download via USB
Expansion Device	Cap tube	Interactive Temperature Graph	Yes
Evaporator Type	Air-cooled fin and tube	Battery Back-up	20 hour display and alarm battery back-up for i.C ³ ® touchscreen (rechargeable 12V, lead acid battery)
Defrost Method	Automatic	Dimensions and Construction	
Refrigerant Charge	R134A Non-CFC	Interior (w x h x d)	19.5 x 17.75 x 22 in 496 x 451 x 559 mm
Performance Data		Exterior (w x h x d) (Includes leveling feet)	24 x 32 x 26.75 in 610 x 813 x 680 mm
Uniformity	+/-1.0°C at 4°C Setpoint	Overall Exterior (w x h x d) (Includes door handle and leveling feet)	24 x 32 x 28 in 610 x 813 x 712 mm
Warranty		ADA Compatible	Fits under 34" (864 mm) ADA compliant counters
Rel.i™ Plus (US/Canada)	7 years compressor, 2 years parts, 1 year labor	Insulation	Sustainable, U.S. EPA and SNAP approved Ecomate®
Outside US/Canada	Consult your local representative	Exterior / Interior Finish	Bacteria-resistant powder coating
		Doors	1, Solid
		Access Ports	Rear access port
		Interior Storage / Capacity	2 perforated metal shelves 18 x 18 in (458 x 458 mm) 100 lb (46 kg) max capacity / shelf
		Leveling Feet	Standard
		Casters	Optional / swivel locking
		Integrated Access Control	Optional - Electromagnetic lock via i.C ³ User Interface PIN
		Net Weight	187 lb (85 kg)
		Shipping Weight	246 lb (112 kg)
		Clearance Requirements	Minimum of 3" (76mm) behind unit for clearance and feature access.
		Options / Accessories	Chart Recorder, Stacking Kit, Floor & Wall Bracket Kit, Remote Alarms, Interior Light, Glass Door, Casters, Stainless Steel Interior and Exterior, Remote Lock Adapter Kit, Extended Warranty



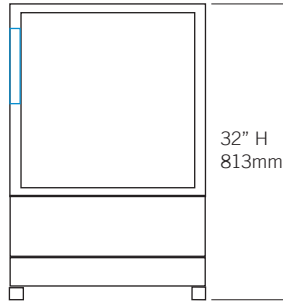

Note: Height with optional casters
 33.5 in (851 mm)

