

Door Handles and Adapters

for Omnicell® FlexLock

Helmer Scientific MedLock Solutions

Special accessories are available for facilities that install the Omnicell FlexLock on Helmer upright refrigerators, undercounter refrigerators, and undercounter freezers.

FlexLock Compatible Door Handles for Upright Models

The door handle features a cutout and spacer which enable the FlexLock to wrap around the front and side of the door.

- Available for single door upright refrigerators
11, 13, 20, and 25 cu ft (326, 377, 572, and 714 L)
- Compatible with glass and solid doors
- Factory or field installation


FlexLock Adapter Kit Undercounter Models

The adapter kit includes an extension bar that mounts to the door of the unit, allowing the FlexLock to be properly installed.

- Adapter kits available for undercounter refrigerators and freezers (5 cu ft / 150 L)
- Compatible with glass and solid doors; powder coated and stainless steel exteriors
- Factory installed or field installation kits available



Cut-out of Door Handle for Upright Refrigerators

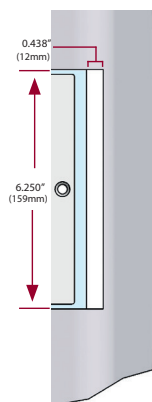

Upright Omnicell®
FlexLock Compatible
Door Handle



Undercounter Adapter
(Shown with Omnicell®
FlexLock)

GX Solutions i.Series® and Horizon Series™ Models

Factory Installed Options

Omnicell FlexLock Adapter Kit, Factory installed, 5 cu ft (150 L)
(for UC Refrigerators and Freezers)

Omnicell FlexLock Compatible Door Handle, Factory installed, 20 and 25 cu ft
(572 and 714 L) (for Refrigerators and Freezers)

Field Installed Options

Omnicell FlexLock Adapter Kit, Field installation, 5 cu ft (150 L)
(for UC Refrigerators and Freezers)

Omnicell FlexLock Compatible Door Handle, Field installation, 11 and 13 cu ft
(326 and 377 L) (for Refrigerators)

Omnicell FlexLock Compatible Door Handle, Field installation, 20 and 25 cu ft
(572 and 714 L) (for Refrigerators and Freezers)

Note:

Adapter kits and compatible door handles do not include the medication dispensing lock.

Flexlock compatible door handles and lock adapter kits also available for legacy refrigerator and freezer models. Contact your sales rep for more information.

*Confirm dimensions of your lock prior to placing order